

## Monolite Trays



**BEST QUALITY**  
  
MADE IN GERMANY

### MONOLITE TRAYS

#### Safe transport and storage solutions

The Monolite trays are transport and storage systems for endoprosthetics, traumatology, dental and other medical areas and are certified as a class I medical device.

The trays are made of one piece, without any material overlapping, from high quality stainless steel. That ensures plenty of thermal capacity during reprocessing, so residual humidity after the cleaning process is not an issue. Their high thermal and chemical resistance and the use of the transport lid as a cover, make transport and storage particularly easy.

The tray system is characterised by its light weight and can be stacked as desired. The individual positioning of the retaining elements in the whole matrix ensure a smooth workflow during any surgical procedure.




THE BENEFITS AT A GLANCE:

- suitable for a variety of container systems
- trays can be delivered empty or complete with fixtures
- individual layout with a variety of fixtures
- individual layout on the bottom of the tray

Individual marking of the tray and tray lid with your company logo is possible.

FOLLOWING DIMENSIONS ARE AVAILABLE:


- |                |                      |                     |                      |
|----------------|----------------------|---------------------|----------------------|
| ▪ Monolite 60  | 510 x 255 x 60,5 mm  | ▪ Mini Monolite 60  | 251 x 255 x 60,5 mm  |
| ▪ Monolite 80  | 510 x 255 x 80,5 mm  | ▪ Mini Monolite 80  | 251 x 255 x 80,5 mm  |
| ▪ Monolite 110 | 510 x 255 x 110,5 mm | ▪ Mini Monolite 110 | 251 x 255 x 110,5 mm |

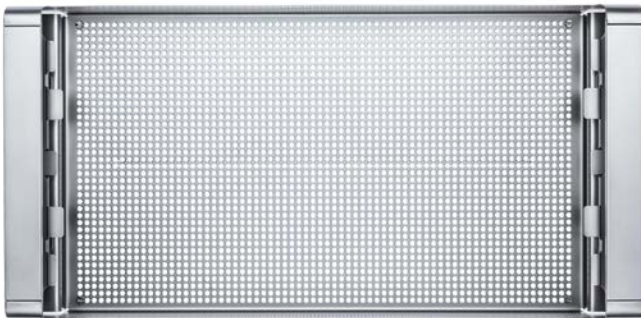

Monolite Tray 110



Mini Monolite Tray 60

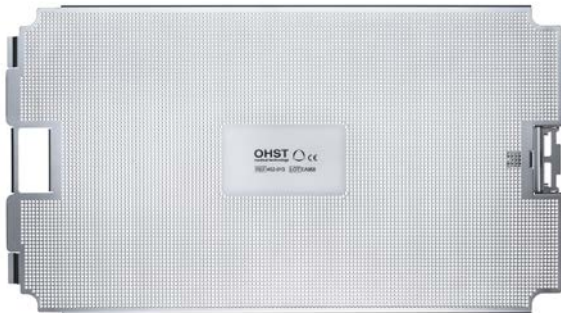

Lid for Monolite Tray 110



Lid for Mini Monolite Tray 60

